




DESS

> DENTAL SMART SOLUTIONS

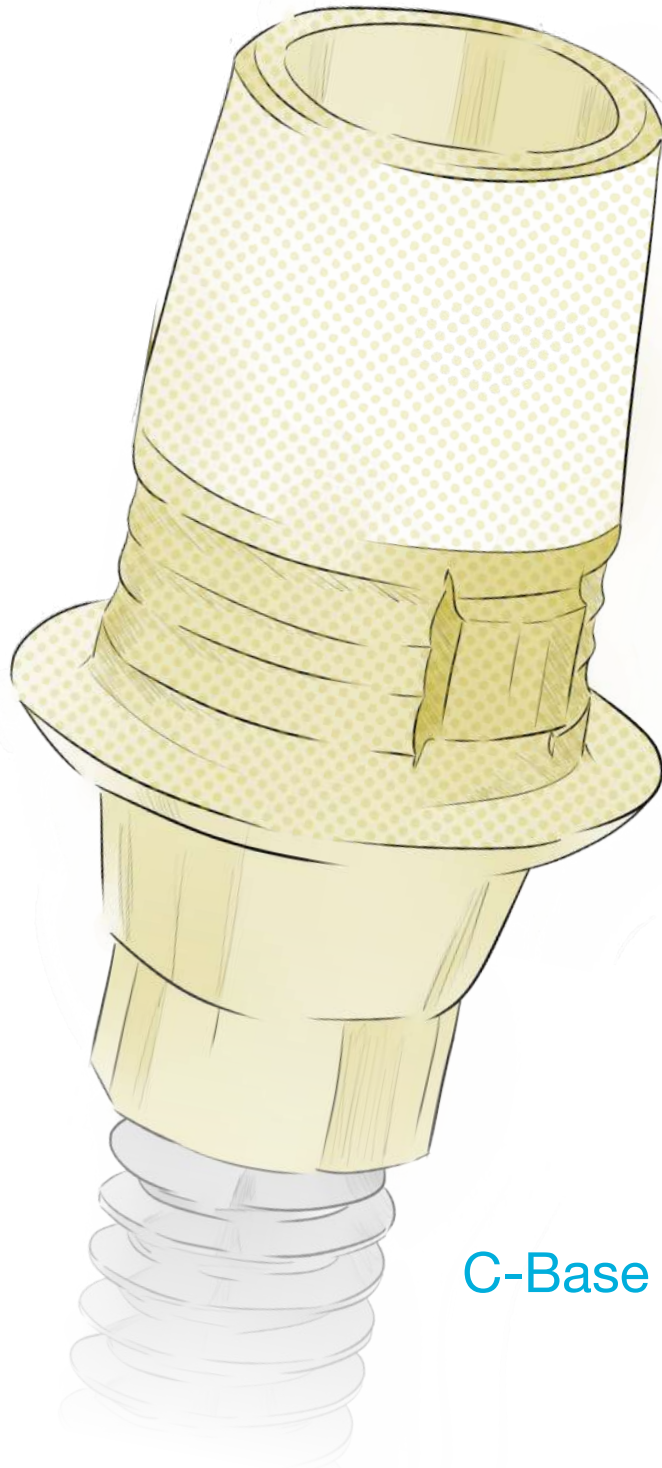

C-Base 2021/1



The Cerec® Sirona® system is undoubtedly one of the most complete and versatile solutions in the market.

Few solutions offer the agility the Cerec® System does, where you can go from scan to finish the piece in one session.

We create a compelling alternative that, based on the “Pure Switch” concept, would offer an improved solution to be used under the Sirona digital workflow with Sirona® libraries or available in DESS® libraries for 3Shape, Dental Wings & Exocad.

For this reason, we included slight modifications that enhance the overall properties of the product, without disrupting the usage.

These modifications can be summarized in:






- Gold anodized colour
- Multiple gingival heights
- SelectGrip® enhanced retention surface treatment
- Screw included, same as implant brand manufacturer
- Non-engaging versions

Manufactured in Titanium Grade 5 ELI, and with 15 implant compatibilities and more than 180 references, the DESS® C-Base® range is the most comprehensive offer of compatible Cerec® Ti-bases in the market.








Digital Workflow

Cerec® / Inlab Systems:

-  **1 Scan either intraoral or on the model**
Use either the Sirona scan post in the C-base or the “S/L” scan bodies
-  **2 Select the Cerec® / Inlab software reference**
To match the implant connection & platform compatibility
-  **3 Follow the Sirona® CAD Software instructions**
Design your restoration as per your needs
-  **4 Use the “S/L” Blocks to mill your abutment**
-  **5 Cement the zirconia mesostructure to the C-Base®**



ExoCad or 3Shape:

-  **1 Scan either intraoral or on the model**
Use the DESS® Desktop or intraoral scan bodies
-  **2 Select the C-Base® reference in ExoCad or 3Shape**
Search for the implant/platform compatibility
-  **3 Design your restoration**
-  **4 Mill your abutment or send to mill**
-  **5 Cement the zirconia mesostructure to the C-Base®**

Reengineered to perfection

Improved aesthetic finish

Gold Anodized colour.

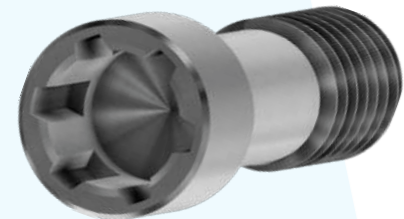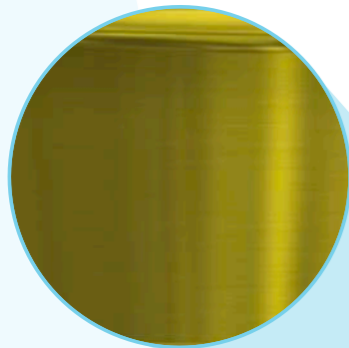


DESS® C-Base includes final screw

Identical to the Implant brand screw.



Engaging or non-engaging

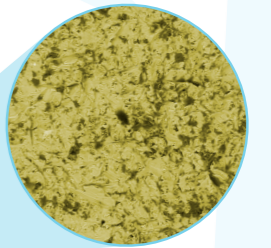

*Microscope view

SelectGrip® surface treatment enhancement, improves overall retention up to **500%**.



Multiple Gingival Heights

Up to 3 different alternatives.



NobelActive® & NobelReplace® CC Active Hex



		NP PURPLE 3.5	RP YELLOW 4.3-5.0	WP
0,8mm	●	58.041	58.042	
3,0mm	●	58.041/3	58.042/3	
0,8mm	●	59.041	59.042	
3,0mm	●	59.041/3	59.042/3	
Replacement Screw		19.027	19.030	

NobelReplace Select™ Tri-lobe



		PURPLE 3.5	YELLOW 4.3	BLUE 5.0	6.0
0,4mm	●	58.004		58.006	58.097
0,5mm	●		58.005		
1,5mm	●	58.004/2	58.005/2	58.006/2	
3,0mm	●	58.004/3	58.005/3	58.006/3	
0,4mm	◀	59.004		59.006	59.097
0,5mm	◀		59.005		
1,5mm	◀	59.004/2	59.005/2	59.006/2	
3,0mm	◀	59.004/3	59.005/3	59.006/3	
Replacement Screw		19.004	19.005	19.005	19.005

Nobel Brånemark® External Hex Universal



		3.5	4.1
1,5mm	●	58.001/2	58.002/2
3,0mm	●	58.001/3	58.002/3
0,6mm	●	59.001	59.002
1,5mm	●	59.001/2	59.002/2
3,0mm	●	59.001/3	59.002/3
Replacement Screw		19.027	19.028

3i Osseotite® External HEX USA



		PURPLE 3.4	BLUE 4.1	YELLOW 5.0
1,5mm	●	58.011/2	58.012/2	58.013/2
3,0mm	●	58.011/3	58.012/3	58.013/3
0,6mm	●	59.011	59.012	59.013
1,5mm	●	59.011/2	59.012/2	59.013/2
3,0mm	●	59.011/3	59.012/3	59.013/3
Replacement Screw		19.034	19.034	19.034

3i Certain® Internal HEX "CLICK"



		NP PURPLE 3.4	RP BLUE 4.1	WP YELLOW 5.0
0,4mm	●	58.014	58.015	
0,5mm	●			58.016
1,5mm	●	58.014/2	58.015/2	58.016/2
3,0mm	●	58.014/3	58.015/3	58.016/3
0,4mm	●	59.014	59.015	
0,5mm	●			59.016
1,5mm	●	59.014/2	59.015/2	59.016/2
3,0mm	●	59.014/3	59.015/3	59.016/3
Replacement Screw		19.011	19.011	19.011

Straumann® Soft Tissue Level & synOCTA® Octagon



		RN 4.8	WN 6.5
1,0mm	●	58.009	
1,2mm	●		58.045
1,0mm	●	59.009	
1,2mm	●		59.045
Replacement Screw		19.026	19.026

Straumann® Bone Level Conical BL



		NC 3.3	RC 4.1
0,7mm	●	58.043	58.044
3,0mm	●	58.043/3	58.044/3
0,7mm	●	59.043	59.044
3,0mm	●	59.043/3	59.044/3
Replacement Screw		19.032	19.033

Zimmer Screw-vent®
Internal Hex USA



		NP GREEN 3.5	RP PURPLE 4.5	WP YELLOW 5.7
0,5mm	●	58.017	58.018	58.019
1,5mm	●	58.017/2	58.018/2	58.019/2
3,0mm	●	58.017/3	58.018/3	58.019/3
0,5mm	●	59.017	59.018	59.019
1,5mm	●	59.017/2	59.018/2	59.019/2
3,0mm	●	59.017/3	59.018/3	59.019/3
Replacement Screw		19.012	19.012	19.012

Astra Tech implant system™ EV
Conic EV



		PURPLE 3.6	YELLOW 4.2	BLUE 4.8
1,0mm	●	58.059	58.060	58.061
3,0mm	●	58.059/3	58.060/3	58.061/3
1,0mm	●	59.059	59.060	59.061
3,0mm	●	59.059/3	59.060/3	59.061/3
Replacement Screw		19.059	19.060	19.061

Astra Tech Osseospeed™
Internal Hex Conic



		AQUA 3.5-4.0	LILAC 4.5-5.0
0,8mm	●	58.024	
0,9mm	●		58.025
3,0mm	●	58.024/3	58.025/3
0,8mm	●	59.024	
0,9mm	●		59.025
3,0mm	●	59.024/3	59.025/3
Replacement Screw		19.016	19.017

Dentsply Friadent Xive®
Internal Hex FD



		SILVER 3.4	YELLOW 3.8	BLUE 4.5	RED 5.5
0,3mm	●	58.038	58.039		
0,6mm	●			58.040	58.071
1,5mm	●	58.038/2	58.039/2	58.040/2	58.071/2
3,0mm	●	58.038/3	58.039/3	58.040/3	58.071/3
0,3mm	●	59.038	59.039		
0,6mm	●			59.040	59.071
1,5mm	●	59.038/2	59.039/2	59.040/2	59.071/2
3,0mm	●	59.038/3	59.039/3	59.040/3	59.071/3
Replacement Screw		19.024	19.024	19.024	19.024

Dentsply Ankylos® C/X
Internal ANK



1,4mm	●	58.050
2,5mm	●	58.050/2
1,4mm	○	59.050
2,5mm	○	59.050/2
Replacement Screw		19.046

OSSTEM® TS
Conic OSS



		Mini YELLOW	STD GREEN regular
1,0mm	●	58.072	58.073
3,0mm	●	58.072/3	58.073/3
1,0mm	●	59.072	59.073
3,0mm	●	59.072/3	59.073/3
Replacement Screw		19.072	19.073

CONELOG®
Internal CON



		NP SILVER 3.3	RP YELLOW 3.8/4.3	WP BLUE 5.0
1,0mm	●	58.074	58.075	58.076
2,1mm	●	58.074/2	58.075/2	58.076/2
1,0mm	●	59.074	59.075	59.076
2,1mm	●	59.074/2	59.075/2	59.076/2
Replacement Screw		19.052	19.052	19.053

BioHorizons® Internal
Internal Hex BH



		NP YELLOW 3.5	RP GREEN 4.5	WP BLUE 5.7
0,5mm	●	58.087	58.088	58.089
3,0mm	●	58.087/3	58.088/3	58.089/3
0,5mm	●	59.087	59.088	59.089
3,0mm	●	59.087/3	59.088/3	59.089/3
Replacement Screw		19.043	19.043	19.043

CONNECTION	SOFTWARE	PRODUCT REFERENCES						SCREW	TORQUE	PLATFORM	
		GH1		GH2		GH3					
		○	○	○	○	○	○				
NOBEL ACTIVE® & NOBEL REPLACE® CC	L	NB A 4.5 L	58.041	59.041			58.041/3	59.041/3	19.027		NP
		NB A 5.0 L	58.042	59.042			58.042/3	59.042/3	19.030	35 Ncm	RP
NOBEL REPLACE SELECT™	L	NB RS 3.5 L	58.004	59.004	58.004/2	59.004/2	58.004/3	59.004/3	19.004		NP (3.5)
		NB RS 4.3 L	58.005	59.005	58.005/2	59.005/2	58.005/3	59.005/3	19.005	35 Ncm	RP (4.3)
		NB RS 5.0 L	58.006	59.006	58.006/2	59.006/2	58.006/3	59.006/3	19.005		WP (5.0)
		NB RS 6.0 L	58.097	59.097	58.097/2	59.097/2	58.097/3	59.097/3	19.005		6.0
NOBEL BRÅNEMARK®	L	NB B 3.4 L		59.001	58.001/2	59.001/2	58.001/3	59.001/3	19.027	35 Ncm	NP (3.4)
		NB B 4.1 L		59.002	58.002/2	59.002/2	58.002/3	59.002/3	19.028		RP (4.1)
3i OSSEOTITE™	L	B 0 3.4 L		59.011	58.011/2	59.011/2	58.011/3	59.011/3	19.034		3.4
		B 0 4.1 L		59.012	58.012/2	59.012/2	58.012/3	59.012/3	19.034	20 Ncm	4.1
		B 0 5.0 L		59.013	58.013/2	59.013/2	58.013/3	59.013/3	19.034		5.0
3i CERTAIN®	S	B C 3.4 S	58.014	59.014	58.014/2	59.014/2	58.014/3	59.014/3	19.011		3.4
		B C 4.1 L	58.015	59.015	58.015/2	59.015/2	58.015/3	59.015/3	19.011	20 Ncm	4.1
		B C 5.0 L	58.016	59.016	58.016/2	59.016/2	58.016/3	59.016/3	19.011		5.0
STRAUMANN® SOFT TISSUE LEVEL & synOcta®		S SO 4.8 L	58.009	59.009					19.026		RN
		S SO 6.5 L	58.045	59.045					19.026		WN
STRAUMANN® BONE LEVEL	L	S BL 3.3 L	58.043	59.043			58.043/3	59.043/3	19.032	35 Ncm	NC
		S BL 4.1 L	58.044	59.044			58.044/3	59.044/3	19.033		RC
ZIMMER® SCREW-VENT®	L	Z TSV 3.5 L	58.017	59.017	58.017/2	59.017/2	58.017/3	59.017/3	19.042		3.5
		Z TSV 4.5 L	58.018	59.018	58.018/2	59.018/2	58.018/3	59.018/3	19.042	30 Ncm	4.5
		Z TSV 5.7 L	58.019	59.019	58.019/2	59.019/2	58.019/3	59.019/3	19.042		5.7
ASTRA TECH IMPLANT SYSTEM™ EV	S	AT EV 3.6 GH1 S	58.059	59.059			58.059/3	59.059/3	19.059		3.6
		AT EV 4.2 GH1 L	58.060	59.060			58.060/3	59.060/3	19.060	25 Ncm	4.2
		AT EV 4.8 GH1 L	58.061	59.061			58.061/3	59.061/3	19.061		4.8
ASTRA TECH OSSEOSPEED™	L	AT OS 3.5/4.0 L	58.024	59.024			58.024/3	59.024/3	19.016	20 Ncm	3.5/4.0
		AT OS 4.5/5.0 L	58.025	59.025			58.025/3	59.025/3	19.017	25 Ncm	4.5/5.0
DENTSPLY FRIADENT XIVE®	S	FX 3.4 S	58.038	59.038	58.038/2	59.038/2	58.038/3	59.038/3	19.024		3.4
		FX 3.8 S	58.039	59.039	58.039/2	59.039/2	58.039/3	59.039/3	19.024	24 Ncm	3.8
	L	FX 4.5 L	58.040	59.040	58.040/2	59.040/2	58.040/3	59.040/3	19.024		4.5
		FX 5.5 L	58.071	59.071	58.071/2	59.071/2	58.071/3	59.071/3	19.024		5.5
DENTSPLY ANKYLOS™ C/X	S	ANKN C/ GH1 S	58.050	59.050					19.046		
		ANKN C/ GH2 S			58.050/2	59.050/2			19.046	15 Ncm	C/X
OSSTEM® TS	L	O TS 3.5 L	58.072	59.072			58.072/3	59.072/3	19.072	20 Ncm	MINI
		O TS 4.0 L	58.073	59.073			58.073/3	59.073/3	19.0737	30 Ncm	STANDARD
CONELOG®	S	CONELOG 3.3/x	58.074	59.074	58.074/2	59.074/2			19.052		3.3
		CONELOG 3.8/x	58.075	59.075	58.075/2	59.075/2			19.052	20 Ncm	3.8/4.3
		CONELOG 5.0/x	58.076	59.076	58.076/2	59.076/2			19.053		5.0
BIOHORIZONS™ INTERNAL	L	BH 3.5 L	58.087	59.087			58.087/3	59.087/3	19.043		3.5
		BH 4.5 L	58.088	59.088			58.088/3	59.088/3	19.043	30 Ncm	4.5
		BH 5.7 L	58.089	59.089			58.089/3	59.089/3	19.043		5.7

All DESS products are manufactured according to ISO 9001 and ISO 13485 with CE marking for all products Class IIb

All DESS® AURUMBase® screws and screwdrivers are based on the standard ISO 10664

3i® and 3i Certain®, Miniplant®, MicroMiniplant™ and Osseofite® are registered trademarks of 3i Implant Innovations, USA

ASTRA TECH™, OsseoSpeed™, ANKYLOS® C/X, XIVE®, FRIADENT®, ASTRATECH implant system™ EV, Cerec®, Sirona® and UniAbutment® are registered trademarks of DENTSPLY SIRONA

BioHorizons® is a registered trademark of BioHorizons, INC.

Brånemark System®, NobelReplace® and Replace Select™, NobelActive™, NobelReplace® CC, NobelSpeedy™ and Multi-unit® are registered trademarks of Nobel Biocare AB, Sweden

ITI®, Straumann® and synOcta® are registered trademarks of Straumann Holding AG, Switzerland

MIS® is a registered trademark of MIS Implant Technologies Ltd.

Zimmer Screw-Vent® and SwissPlus® are registered trademarks of Zimmer Dental INC, USA

Dyna® is a registered trademark of Dyna Dental Engineering BV

Locator® is a registered trademark of ZEST IP Holdings

Camlog® and Conelog® are registered trademarks of CAMLOG Biotechnologies AG

Megagen® is a registered trademark of MEGAGEN IMPLANT CO. LTD

TORX® is a trademark of ACUMENT INTELLECTUAL PROPERTIES, LLC

Bego Semados® is a registered trademark of BEGO GmbH & Co. KG

Anthogyr® is a registered trademark of Anthogyr SAS

Medentis ICX® is a registered trademark of medentis medical GmbH

Osstem® is a registered trademark of OSSTEM IMPLANT CO. LTD

Dentium® is a trademark of Dentium Co. Ltd.

GOBALD® is a trademark of Global D SASU.

Biotech Kontakt® is a trademark Biotech Dental SAS.

BTI® is a trademark BTI Biotechnology Institute.

PHIBO® is a trademark PHIBO Dental Solutions S.L.

DESS products can be used only by licensed dentists, physicians and labs.
Please refer ever to the instructions of use available on the website
www.dess-abutments.com
or supplied with the product.

Some products may not be available in your country. Some of the products are cleared for US and Canada.
Please contact your local distributor for more details.



